

On-Call

Express Voice

Blood Glucose **Monitoring System**

Part of the On Call Express Family



**VOICE ENABLED
FOR EASY USE**

Speaks in English
and Spanish

Convenience, Accuracy and Reliability

- ▶ Convenient, affordable and audible Blood Glucose Monitoring System
- ▶ Reliable testing with clinically proven accuracy using On Call Express test strips
- ▶ Large display for maximum visibility
- ▶ 24 hour On Call Customer Support 1-800-838-9502

On Call® Express Voice – Reliable, Easy and Convenient

On Call Express Voice is a convenient, affordable and fully audible Blood Glucose Monitoring System.

Its tactile features are designed for easy navigation and audible instructions to assist the user.

The small blood sample size makes it easy for a patient to utilize alternate site testing. Furthermore, it is based on the same trusted advanced biosensor technology as all other On Call Express Blood Glucose Monitoring Systems.

On Call Express Voice is the perfect solution to glucose monitoring system when convenience matters.



On Call Express Voice

- English and Spanish audible instructions assist the user from set up to testing.
- Tactile features are designed for easy navigation.
- The high contrast strip port is easy to see and shaped for effortless insertion.
- Large display with bold numbers for maximum visibility

Reliable


- ✓ Clinically Proven Accuracy
- 🗣️ Speaks in English and Spanish
- ⚠️ Hypo, Hyper and Ketone Alerts
- 🕒 5 Daily Testing Reminders
- 🍴 Meal Markers

Fast and Easy

- 🕒 4 Second Testing
- 💧 Only 0.4 μ L Blood Sample
- 📍 Alternate Site Testing
(Palm, Forearm)
- 🚫 No Coding Required

Advanced Biosensor Technology for Fast and Accurate Glucose Testing

- ▶ Auto calibration circuitry inside every strip
- ▶ Accurate results in just 4 seconds with only a 0.4 µL blood sample
- ▶ Less painful testing from palm and forearm

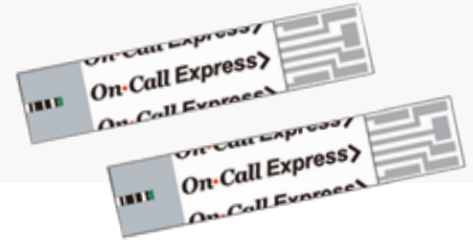

Voice Meter Comparison Chart

					
Company	ACON Laboratories	Apex Bio	OK Biotech	BroadMaster Biotech	i-SENS
System	On Call® Express Voice	Omnis Embrace®	Prodigy Autocode®	Advocate® Redi-Code+Speaking	CareSens® N voice
Enzyme	Glucose Oxidase	Glucose Oxidase	Glucose Oxidase	Glucose Oxidase	Glucose Oxidase
Blood Sample	Capillary	Capillary	Capillary	Capillary	Capillary
Sample Volume	0.4 µL	0.6 µL	0.7 µL	1.1 µL	0.5 µL
Test Time	4 sec.	6 sec.	7 sec.	5 sec.	5 sec.
Meas. Range	20 - 600 (mg/dl)	20 - 600 (mg/dl)	20 - 600 (mg/dl)	20 - 600 (mg/dl)	20 - 600mg/dl
Hematocrit Range	25 - 60%	30 - 55%	30 - 55%	20 - 60%	25 - 60%
Strip Ejector	Yes	No	No	No	Yes
Languages	English, Spanish	English, Spanish	English, Spanish, French, Arabic	English, Spanish	English, Spanish
Memory	500 tests	300 tests	450 tests	400 tests	500 tests
Size (LxWxH)	3.9" x 2.4" x 0.9"	3.9" x 2.1" x 0.8"	3.8" x 2.1" x 0.8"	3.7" x 2.5" x 1.1"	3.9" x 2.2" x 0.7"
Weight	3.2 oz	2.3 oz	2.8 oz	1.9 oz	2.5 oz

Specifications

FEATURE	SPECIFICATION
Measurement Range	20 – 600 mg/dL
Result Calibration	Plasma-equivalent
Sample	Fresh capillary whole blood
Minimum Sample Size	0.4 µL
Test Time	4 seconds
On/Off Source	2 AAA batteries
Battery Life	1,000 tests
Glucose Units of Measure	The meter is pre-set to milligrams per deciliter (mg/dL)
Memory	Up to 500 records with time and date
Automatic shutoff	2 minutes after last action or 5 sec after strip removal
Meter Size (LxWxH)	3.9" x 2.4" x 0.9"
Display Size (LxW)	1.5" x 1.7"
Operating Temperature	41 – 113 °F
Operating Relative Humidity	10 – 90% (non-condensing)
Hematocrit Range	25 – 60%
Data Transfer	USB

Ordering Information

PRODUCT NAME	CATALOG NO.	NDC	CONTENTS
On Call® Express Voice Blood Glucose Monitoring System	757101	82607-0571-01	1 On Call Express Voice Meter, 10 On Call Express test strips, Control Solution 1, 10 Lancets (30 G), 1 On Call Stealth Lancing Device w/ clear cap, 1 Carrying Case, User's Manual in English and Spanish, Quick Reference Guide (En/ Es), Quick Start Guide (En/Es), 1 Warranty Card (En/Es), 1 Logbook, 2 AAA Batteries pre-installed
On Call Express Voice Blood Glucose Meter	757102	82607-0571-02	1 On Call Express Voice Meter, 1 Carrying Case, User's Manual (En/Es), Quick Reference Guide (En and Es), Quick Start Guide (En/Es), 1 Warranty Card (En/Es), 2 AAA Batteries pre-installed
On Call Express Test Strips	755729	82607-0557-29	50 test strips (50/vial), package insert
On Call Express Glucose Control Solution	755732	82607-0557-32	2 mL control solution 1, 2 mL control solution 2, 2 mL control solution 3, package insert
On Call Stealth Lancing Device	754182	82607-0541-82	Lancing device, clear cap for AST, package insert
On Call Sterile Lancets	754105	82607-0541-05	100 sterile 30G lancets (25/bag)



aconlabs.com

ACON Laboratories, Inc.
10125 Mesa Rim Road, San Diego, CA 92121, U.S.A.
Tel: 1.800.838.9502
Fax: 1.858.200.0729
Email: info@aconlabs.com