



Strep A Rapid Test Cassette

A rapid test for the qualitative detection of Strep A antigen in throat swab specimens.



Rapid Testing



Easy to Use



Reliable Results

Flowflex Plus Strep A Rapid Test Cassette

The Flowflex Plus Strep A Rapid Test Cassette is a rapid chromatographic immunoassay for the qualitative detection of Group A Streptococcus antigen from throat swab specimens from symptomatic patients. The test is used to aid in the diagnosis of Group A Streptococcal infection. All negative test results should be confirmed by bacterial culture because negative results do not preclude infection with Group A Streptococcus and should not be used as the sole basis for treatment.

- Throat swab specimens
- Result in 5 minutes
- 24 months shelf life
- Store between 36 to 86°F
- Excellent performance compared to lab-based throat cultures
- CLIA waived
- Includes 2 additional test cassettes for QC testing


Clinical Performance

The Flowflex Plus Strep A Rapid Test Cassette was evaluated against bacterial culture at a CLIA-certified laboratory. The test correctly identified 98.2% of positive specimens and 95.1% of negative specimens.

Materials Provided

- Test Cassettes
- Test Tubes
- Dropper Tips
- Throat Swabs
- Package Insert
- Quick Reference Instructions
- Reagent A (2M Sodium Nitrite)
- Reagent B (0.2M Acetic Acid)
- Positive Control (Heat-inactivated Group A Streptococcus)
- Negative Control (Heat-inactivated Group C Streptococcus)

Test Procedure and Interpretation

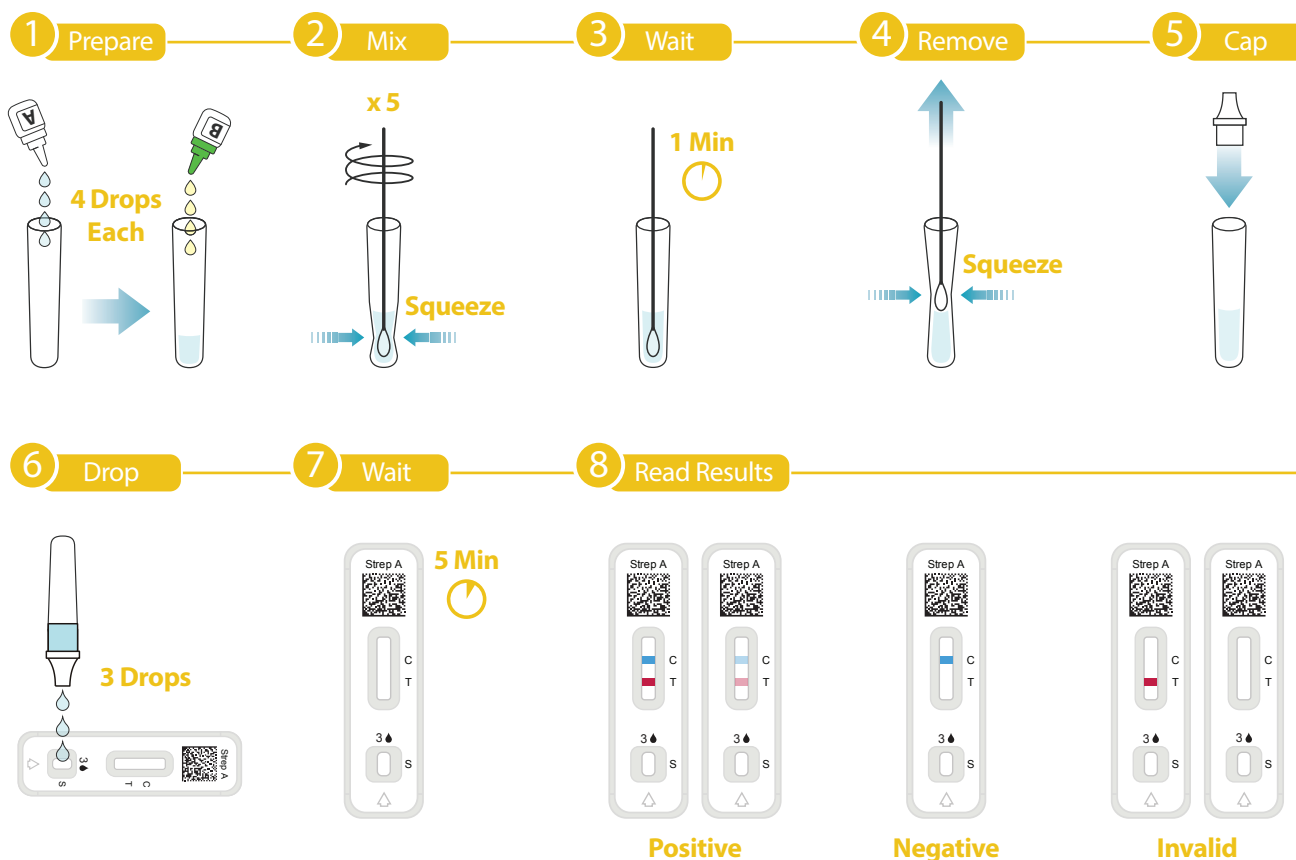

Ordering Information

Product Name	Catalog No.	Format	Specimen	Package
Flowflex Plus Strep A Rapid Test Cassette	L03A-R0115	Cassette	Throat Swabs	25 Tests/Kit

ACON Laboratories, Inc. 9440 Carroll Park San Diego, CA 92121 | Tel: 1.858.875.8000 Fax: 1.858.200.0729 | Email: info@aconlabs.com